

ID-Nr. 30612-001
Model SL.196.060.MH



Cabinet for acids and alkalis SL-CLASSIC
Model SL.196.060.MH, light grey (RAL 7035)

Safe and approved storage of aggressive, non flammable hazardous materials in working areas

Function / construction:

- **Robust construction and longevity:** cabinet body made of powder-coated sheet steel, safety elements assembled outside the storage compartment for increased protection against corrosion
- **Two in one:** two hermetically separated storage compartments allow the isolated storage of acids and bases in one cabinet
- **No unauthorised use:** doors lockable with profile cylinder (integration in an existing locking system possible)
- **Easy alignment:** adjusting aids to compensate for uneven floor
- **ventilation:** integrated metal-free air ducts ready for connection to a technical exhaust system with only one extraction air inlet

| Dimensions | |
|-------------------|------------|
| Width (external) | 597.00 mm |
| Depth (external) | 603.00 mm |
| Height (external) | 1965.00 mm |
| Width (internal) | 485.00 mm |
| Depth (internal) | 570.00 mm |
| Height (internal) | 914.00 mm |

The available internal dimensions can vary depending on the interior equipment package.


| Weight / Load | |
|-----------------------------------|--------------------------|
| Weight without interior equipment | 136 kg |
| Distributed load | 444.00 kg/m ² |
| Maximum load | 300 kg |

| Ventilation | |
|---|---------------------|
| Extraction air | 75 DN |
| Recom. exhaust air quantity 30 times | 8 m ³ /h |
| Differential pressure (30 times air change) | 40 Pa |
| Minimum exhaust air quantity 10 times | 5 m ³ /h |
| Differential pressure (10 times air change) | 40 Pa |

| Colouring / Surface finish | |
|------------------------------|------------|
| Body colour code | RAL 7035 |
| Body colour | light grey |
| Body surface finish | plain |
| Door / Drawer colour code | RAL 7035 |
| Door / Drawer colour | light grey |
| Door / Drawer surface finish | plain |

| General Information | |
|-----------------------|----------|
| Number of doors | 2 |
| Depth with open doors | 1119 mm |
| Customs tariff number | 94032080 |

Front view

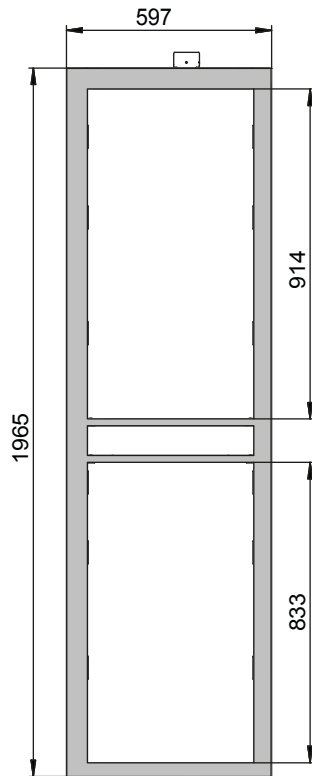

Side view

